

**LECKEY<sup>®</sup>**

# Squiggles Seating System

~~Star~~  
Squiggles



# LECKEY<sup>o</sup>

Stronger together.

We work with individuals, therapists and carers to design products with both a clinical and an emotional focus. Using the latest research and clinical understanding, we create practical solutions which are easily integrated into family life, because for us, life is about going, enjoying, participating and doing.

# It's all about progress.

24 hour  
postural care for  
babies, kids & adults.  
Sleeping, Sitting,  
Standing, Walking,  
Moving, Bathing,  
Toileting.





We are dedicated to ensuring we play our part in helping kids to make the most out of every day by following our key design principles:

- 1 Child-centred & clinically sound
- 2 Easy-to-use & lightweight
- 3 Durable & dependable
- 4 Attractive & desirable
- 5 Recyclable & reusable



  
Squiggles





The Squiggles Seating System for children aged 1-5 years benefits from a number of improved design features and more advanced materials, making it more comfortable and easier to care for.

#### **Practicality**

Removeable machine washable covers, 40°C.

#### **Durability**

Allen key screws supplied as standard with anti-vibration washers to prevent loosening through continued use.

#### **Comfort**

A 30mm one piece seat base cushion ensures comfort, allowing the child to focus on other activities.



**And we offer so much more!**

#### **Head**

Wide, contoured headrest providing greater lateral support.

#### **Trunk**

Contoured shoulder pad increases back height by 75mm.

Easy-to-use tray locking mechanism.

Padded grab rail can be used instead of tray in buggies.

#### **Pelvis**

Pelvic harness attachment points and webbing reinforced to cater for kids with strong extensor patterns.

Width adjustment on hip laterals up to 260mm.

#### **Legs & Feet**

Additional lateral padding at knees, up to 30mm both sides. Helps prevent low tone kids' legs from abducting - accessory only.

#### **Chassis**

Tilt-in-space on Hi-low chassis, giving 30° recline.






With its bright colours and clinical focus, the Squiggles Seating System has always been about fulfilling the postural and emotional needs of young children and their families. The Squiggles Seating System is able to help a broad range of kids to carry out everyday functions at home, in the classroom or even out and about.





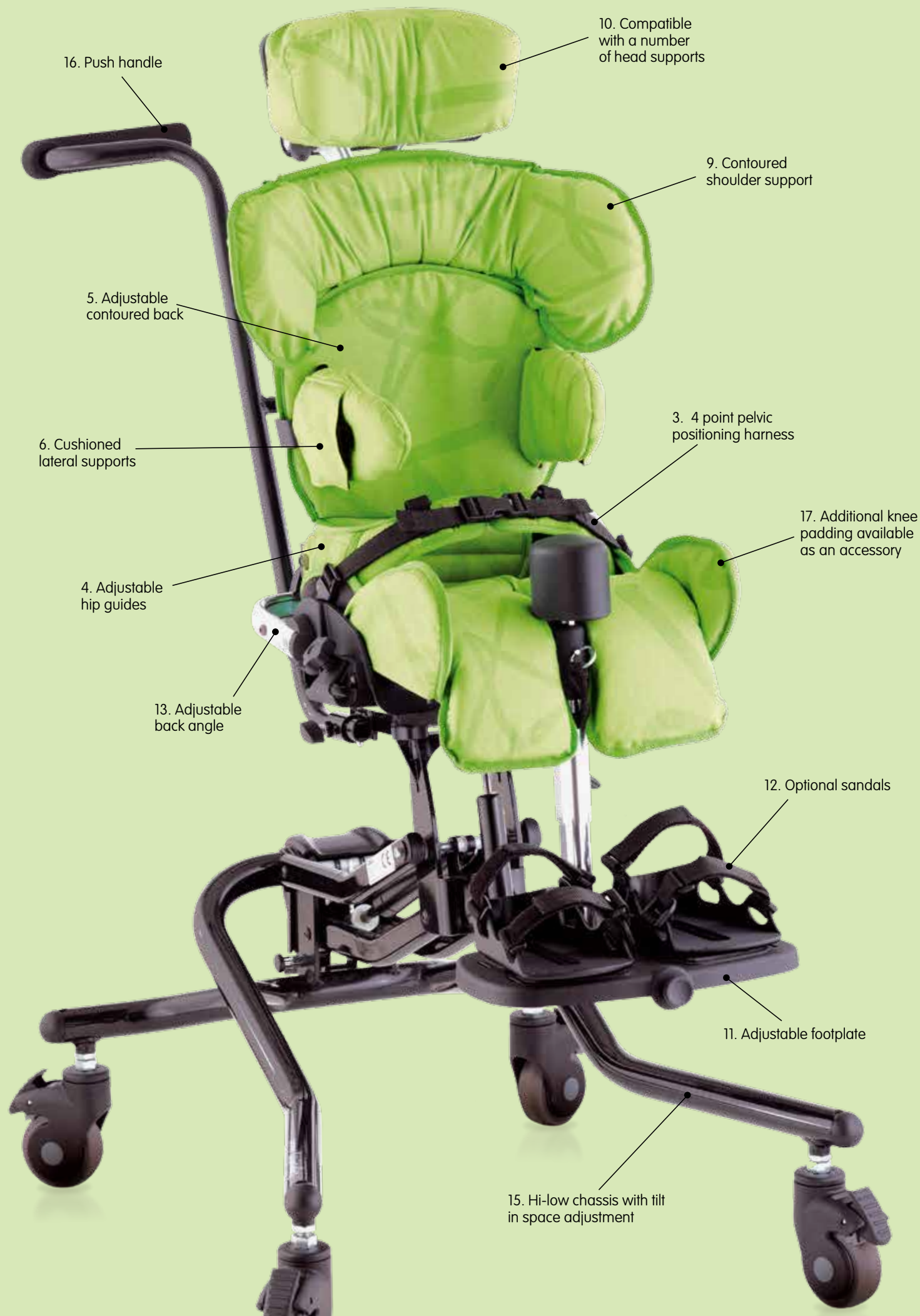

1. Flexible sacral cushion, which can be shaped to the child's lower back, supports the lower spine in the desired position by encouraging a degree of forward tilt in the pelvis and lower trunk extension.
2. Ramped base cushion prevents the child sliding forward.
3. 4 point pelvic harness is adjustable to give a secure, stable base of support and is cushioned for the child's comfort.
4. Hip guides are also adjustable for optimised postural support.
5. Contoured back can be adjusted to the child's height.
6. Cushioned lateral supports provide side support to help maintain a safe, upright position.
7. Chest harness can be attached to provide further support whilst maintaining the child's freedom to enjoy activities.
8. Trunk harness may be used if the child would benefit from side, front and shoulder support.
9. Contoured shoulder support angle adjusts to provide shoulder protraction which helps balance the head and trunk over the pelvis.
10. Headrest is compatible with a variety of head supports: Whitmyer, Otto Bock and other leading brands.
11. Adjustable footplate, which can be tilted to the preferred angle, ensures that the feet are well supported, providing a secure base for maximum upper body function. The footplate flips away to allow for easy transfers.
12. Sandals can be used to position the feet securely.
13. Adjustable back angle mechanism maintains the position of the head and trunk supports as the back angle is changed.
14. Activity tray supports a wide range of activities. A grab rail can also be attached to assist the child if required.
15. Hi-low chassis can adjust from floor to table height allowing the child to enjoy a wide range of activities from circle time to family meals. The chassis has angle adjustment for tilt in space supporting various postures. All adjustments can be made safely with the child in the seat.
16. Push handles allow the chair to be easily moved around the home or classroom.



1. Flexible lumbar and sacral cushion



2. Ramped base cushion



7. Optional chest harness



8. Trunk harness



14. Activity tray









Established in 1983, Leckey is a globally recognised pioneer in the research and development of products that help adults and children with disabilities to go, do, enjoy and participate in everyday activities throughout the day and night.

We take a highly clinical approach to product design and development. Through in-depth clinical research studies with leading universities, and extensive trials with occupational therapists, physiotherapists, users and their families, we continue to develop posturally supportive, family friendly products for all day care, at every stage of life.

Through early intervention, childhood and adulthood Leckey's experienced team of designers, therapists and bioengineers work together to develop products that meet the clinical needs of the healthcare professionals and the social needs of the user.

To achieve this, we work with the healthcare professionals, the individuals and carers who use our products everyday. With their help, we create the dependable, durable, proven and high performance products that we are known for worldwide.

## Free Leckey Servicing

Our free service and support includes:

1. Free product training
2. Free product set-up
3. Free product assessment
4. Free product re-assessment
5. Free maintenance within 2 year warranty

For service policies on all products outside warranty, please contact Leckey's Customer Service department. James Leckey Design Ltd as manufacturer with sole responsibility declares that all products conform to 93/42/EEC guidelines and EN12182 technical aids for disabled persons general requirements and test methods. Order forms and spare parts lists to extend the service life of the product and allow reissue are available on request or online at [www.leckey.com](http://www.leckey.com).





24 hour postural  
care for babies,  
kids & adults.  
Sleeping, Sitting,  
Standing, Walking,  
Moving, Bathing,  
Toileting.